Smart Internz

Pneumonia Detection using X-rays using Watson Studio

Index

Topics		Page No
1.	Introduction	1
	1.1. Overview	
	1.2. Purpose	
2.	Literature Survey	2
	2.1. Existing Problem	
	2.2. Proposed solution	
3.	Theoretical Analysis	3
	3.1. Block Diagram	
	3.2. Hardware / Software Designing	
11.	Experimental Investigations Flowchart Result Advantages & Disadvantages Applications Conclusion Future Scope Bibliography Appendix	4 5 6 7 8 9 9 9
	12.1. Source Code	
	12.2. UI Output Screenshot	

1. <u>Introduction</u>:

1.1. Overview

Pneumonia is an infection that inflames the air sacs in one or both lungs. The air sacs may fill the pus causing cough with pus, fever, chills and difficulty in breathing. Patients with pneumonia may need to be hospitalized or admitted to the ICU which is more than a million and over 50,000 deaths each year. It is the world's leading cause of death for children under the age of 5. In US, it is the most common cause of hospitalization for children and adults.

1.2. Purpose

Chest x-rays have been considered as the best tool to detect any form of Pneumonia. Studies have shown that even experienced radiologists have a hard time to correctly identify whether something on the x-rays is an infiltrate because the pus is denser than air. The lack of doctors and treatment of Pneumonia may also become more difficult. So, an accurate and fast diagnosis of Pneumonia is urgently needed, especially for people in poverty. So, the main purpose is to build a model to automatically identify whether a patient is suffering from Pneumonia or not by checking x-ray images.

2. <u>Literature Survey</u>:

2.1. Existing Problem


Lack of doctors or inaccurate predictions of Pneumonia may make the situations worse and may cause deaths especially for poor and uneducated people as they don't know what to do and what doctors say when they are attacked with Pneumonia. Especially, in these tough situations like COVID-19, doctors have no time to check x-rays and detect the Pneumonia as they are busy in serving COVID-19 patients. In some cases, even though a doctor checks an x-ray of a patient who is mildly attacked by the Pneumonia but considered it as non-infected, his situation may become serious in the coming days.

2.2. Proposed Solution

Keeping the problems and difficulties of Pneumonia patients in mind, a new solution is proposed. An app has been built so that the users just need to upload x-rays and click on predict. Then the prediction will be displayed after clicking on predict button. Users can make use of it in these tough situations or when they are unable to consult a doctor or when they want accurate prediction. It is easy to use and free of cost. Based on that prediction, they can move on to further medications or they may consult doctors virtually as checking x-ray physically is not required.

3. Theoretical Analysis:

3.1. Block Diagram

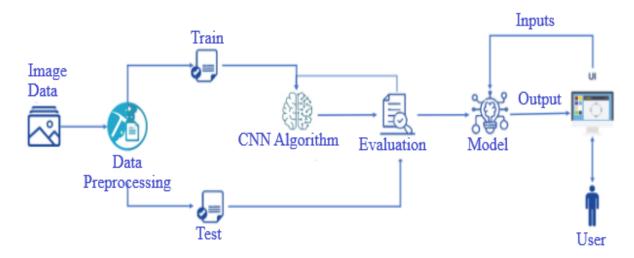

3.2. Hardware / Software Designing

Software Requirements:

- **★** OS Windows XP 7/8/10
- **★** Anaconda
- ★ Jupyter notebook/Pycharm/Spyder
- **★** Python
- ★ Cloud storage space in IBM


Hardware Requirements:

- ★ I3 processor
- ★ RAM 1GB

4. Experimental Investigations:

In this project, the main aim is to predict whether a person is affected with pneumonia or not. So, firstly, image preprocessing is done in which the images are classified into 2 classes – NORMAL, PNEUMONIA. After preprocessing, model is initialized and CNN layers are added to that model and it is fitted and saved. At last, we test the model. After testing the model, an app is built for user interface. Users can access that app and they can upload their x-rays and they themselves can predict whether they are infected or they are alright. So simply and easily, we can predict Pneumonia with an x-ray at our own place and at any time at free of cost.

5. Flowchart:

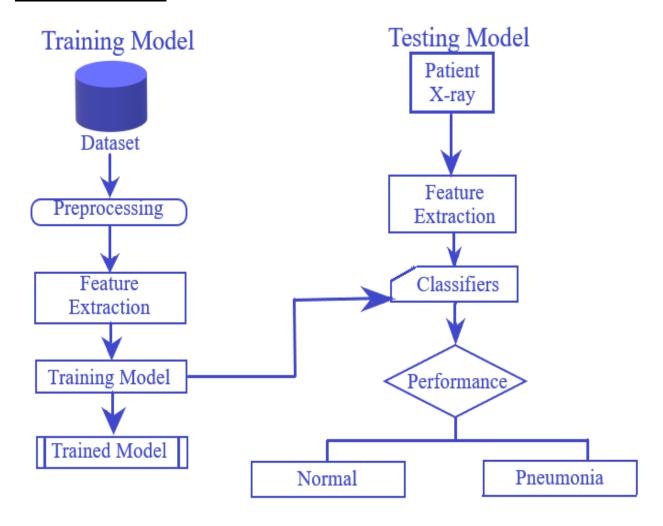

6. Result:

- → The user will be uploading an x-ray image through the UI interface which is created by flask application.
- → The input image is loaded into the model.
- → This model performs the preprocessing and classifies the class name of the uploaded image(0-Normal, 1-Pneumonia) and predicts whether the x-ray is of affected person or not.
- → The output will be sent to the flask app using an http request and when the user clicks on predict, it displays the prediction.

7. Advantages and Disadvantages:

Advantages:

- © Very simple and easy as the user just need to upload the x-ray image.
- © Free of cost. The user can predict pneumonia without consulting doctor.
- © Prediction is 90% accurate.
- © Early prediction of Pneumonia before the situation becomes worse as early prediction may just need medication rather than hospitalization.
- © It is mostly useful for people who is living in remote areas where there may be no doctors.

Disadvantages:

- © For uneducated people and for old people, it may take some time to get used to it.
- © Users must have smart phones and internet to capture the x-ray images and to upload them into the app for prediction.
- © Poor people who do not have at least a smart phone may have no use with this app.

8. Applications:

- ➤ In hospitals there is no need for doctors to predict Pneumonia. The compounders who have basic knowledge may predict and may suggest the patients to consult the doctor for medication.
- ➤ It can be used at residentials houses where people with basic knowledge of mobiles can use it.
- ➤ Today's youth have basic skills about smart phones. So, they can help poor and old people or children who are living in their surroundings to predict the disease.
- ➤ It can be used anywhere and anybody who is having smartphones and internet.

9. <u>Conclusion</u>:

- ✓ Our CNN model can effectively detect pneumonia disease with extremely high accuracy.
- ✓ Radiologists could combine this model and their expertise to make a highly accurate diagnosis of Pneumonia.
- ✓ Mortality rate can be reduced by early prediction of Pneumonia.

10. Future Scope:

In the future, it would be interesting to see approaches in which the weights corresponding to different models can be estimated more efficiently and a model that takes into account the patient's history while making predictions.

11. Bibliography:

- https://towardsdatascience.com/pneumonia-detection-with-keras-and-fastapi-6c10dab657e0
- <u>Diagnostics | Free Full-Text | Efficient Pneumonia Detection in Chest Xray Images Using Deep Transfer Learning | HTML (mdpi.com)</u>
- Chest X-Rays Pneumonia Detection Using Convolutional Neural
 Network | by Max | Analytics Vidhya | May, 2021 | Medium
- o https://www.ncbi.nlm.nih.gov/pmc/articles/PMC7345724/

12. Bibliography:

12.1. Source Code

```
In [1]: from keras.preprocessing.image import ImageDataGenerator, load_img
                   from keras.models import Sequential
                   from keras.layers import Conv2D, MaxPooling2D
                   from keras.layers import Activation, Dropout, Flatten, Dense
                   from keras import backend as K
                   import os
                   import numpy as np
                   import pandas as np
                   import matplotlib.pyplot as plt
                   %matplotlib inline
          In [2]: | train = 'E:/project/dataset/chest_xray/train'
                   test = 'E:/project/dataset/chest_xray/test'
                   val = 'E:/project/dataset/chest_xray/val'
          In [3]: img_width,img_height= 150,150
                   input_shape = (img_width,img_height,3)
In [4]: model = Sequential()
         # The number of filters are 32 and the kernal_size is (3,3) model.add(Conv2D(32, (3, 3), input_shape=input_shape)) model.add(Activation('relu'))
         model.add(MaxPooling2D(pool\_size=(2, 2)))
```

```
In [4]: model = Sequential()
# The number of filters are 32 and the kernal_size is (3,3)
model.add(Conv2D(32, (3, 3), input_shape=input_shape))
model.add(Activation('relu'))
model.add(MaxPooling2D(pool_size=(2, 2)))

model.add(Conv2D(32, (3, 3)))
model.add(MaxPooling2D(pool_size=(2, 2)))

model.add(Conv2D(64, (3, 3)))
model.add(Conv2D(64, (3, 3)))
model.add(MaxPooling2D(pool_size=(2, 2)))

model.add(MaxPooling2D(pool_size=(2, 2)))

model.add(Conv2D(128, (3, 3)))
model.add(Activation('relu'))
model.add(MaxPooling2D(pool_size=(2, 2)))

model.add(MaxPooling2D(pool_size=(2, 2)))

model.add(Flatten())
model.add(Propout(0.5))
model.add(Dropout(0.5))
model.add(Dropout(0.5))
model.add(Dropout(0.5))
model.add(Convation('relu'))
```

```
In [5]: model.compile(loss='binary_crossentropy',
                 optimizer='rmsprop
                 metrics=['accuracy'])
 In [6]: batch_size = 16
         train_datagen = ImageDataGenerator(rescale=1. / 255,shear_range=0.2,zoom_range=0.2,horizontal_flip=True)
        test_datagen = ImageDataGenerator(rescale=1. / 255)
 In [7]: train_generator = train_datagen.flow_from_directory(
           train,
           target_size=(img_width, img_height),
          batch_size=16,
           class_mode='binary')
        test_generator = test_datagen.flow_from_directory(
           target_size=(img_width, img_height),
          batch size=16.
          class_mode='binary')
        validation_generator = test_datagen.flow_from_directory(
           target_size=(img_width, img_height),
          batch_size=16,
          class_mode='binary')
        Found 5216 images belonging to 2 classes.
        Found 624 images belonging to 2 classes.
        Found 16 images belonging to 2 classes.
In [8]: model.fit_generator(
          train_generator,
         steps_per_epoch=5217 // 16,
epochs = 15,
          validation_data=validation_generator,
         validation_steps=17 // 16)
       C:\Users\Home\anaconda3\lib\site-packages\tensorflow\python\keras\engine\training.py:1844: UserWarning: `Model.fit_generator` is deprecate
       d and will be removed in a future version. Please use `Model.fit`, which supports generators. warnings.warn('`Model.fit_generator` is deprecated and '
       Epoch 1/15
       326/326 [==
Epoch 2/15
                            ============] - 1068s 3s/step - loss: 0.6060 - accuracy: 0.7155 - val_loss: 0.6731 - val_accuracy: 0.5625
       326/326 [=====================] - 812s 2s/step - loss: 0.3610 - accuracy: 0.8449 - val_loss: 0.7464 - val_accuracy: 0.8125
       Enoch 3/15
                             ========] - 734s 2s/step - loss: 0.2939 - accuracy: 0.8877 - val_loss: 0.5322 - val_accuracy: 0.8750
       326/326 [====
       Epoch 4/15
                              326/326 [==
       Epoch 5/15
       326/326 [===
                            ========] - 574s 2s/step - loss: 0.2771 - accuracy: 0.9066 - val_loss: 1.2444 - val_accuracy: 0.6250
       Enoch 6/15
                            ==========] - 525s 2s/step - loss: 0.2633 - accuracy: 0.9061 - val_loss: 0.8852 - val_accuracy: 0.7500
       326/326 [====
       Epoch 7/15
                            :==========] - 499s 2s/step - loss: 0.2618 - accuracy: 0.8993 - val_loss: 1.3730 - val_accuracy: 0.6250
       326/326 [===
       Epoch 8/15
       326/326 [==
                               :========] - 481s 1s/step - loss: 0.2404 - accuracy: 0.9148 - val_loss: 1.2890 - val_accuracy: 0.6250
       Epoch 9/15
                           :===========] - 474s 1s/step - loss: 0.3186 - accuracy: 0.9194 - val_loss: 1.2960 - val_accuracy: 0.6250
```

326/326 [=============] - 496s 2s/step - loss: 0.2397 - accuracy: 0.9133 - val_loss: 0.6559 - val_accuracy: 0.7500

326/326 [=== Epoch 10/15

```
In [9]: test_acc = model.evaluate_generator(test_generator,624/16)
print("\nAccuracy:"+" %.2f%%" % (test_acc[1]*100))

C:\Users\Home\anaconda\lib\site-packages\Home\python\keras\engine\training.py:1877: UserWarning: `Model.evaluate_generator` is deprecated and will be removed in a future version. Please use `Model.evaluate`, which supports generators.

warnings.warn('`Model.evaluate_generator` is deprecated and '

Accuracy: 90.38%

In [10]: model.save('mymodel.h5')


In [11]: img_n = load_img('E:/project/dataset/chest_xray/train/NORMAL/IM-0234-0001.jpeg')
plt.imshow(img_n)
plt.ititle("Normal")
plt.show()
```

```
https://www.mg_p)
plt.mshow(img_p)
plt.title("Pneumonia")
plt.show()

Normal
```

```
img_p = load_img('E:/project/dataset/chest_xray/train/PNEUMONIA/person7_bacteria_24.jpeg')
plt.ititle("Pneumonia")
plt.show()
```

 $img_p = load_img('E:/project/dataset/chest_xray/train/PNEUMONIA/person7_bacteria_24.jpeg')$

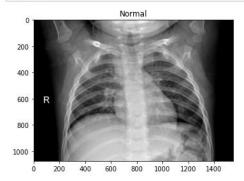




```
In [17]: from tensorflow.python.keras.models import load_model
         from keras.preprocessing import image
         import numpy as np
         import cv2
         model = load_model('mymodel.h5')
model.compile(loss='categorical_crossentropy',
                       optimizer="adam",
                       metrics=['accuracy'])
         from skimage.transform import resize
         def detect(frame):
              #print(type(frame))
           #try:
              img = resize(frame,(150,150))
              img = np.expand_dims(img,axis=0)
              #if(np.max(img)>1):
              #\img = img/255.0
              prediction = model.predict(img)
              print(prediction)
              prediction = model.predict_classes(img)
              print(prediction)
           #except AttributeError:
             # print("shape not found")
In [21]: frame=cv2.imread("E:/project/dataset/chest_xray/test/PNEUMONIA/person109_bacteria_512.jpeg")
         data = detect(frame)
         [[0.9969543]]
[[1]]
In~[22]:~frame=cv2.imread("E:/project/dataset/chest\_xray/test/NORMAL/NORMAL2-IM-0249-0001.jpeg")
         #print(frame)
         data = detect(frame)
         [[0.04447159]]
[[0]]
In [23]: frame=cv2.imread("E:/project/dataset/chest_xray/test/NORMAL/NORMAL2-IM-0290-0001.jpeg")
         #print(frame)
data = detect(frame)
         [[0.01325998]]
         [[0]]
In [24]: frame=cv2.imread("E:/project/dataset/chest_xray/test/PNEUMONIA/person137_bacteria_655.jpeg")
         #print(frame)
         data = detect(frame)
         [[0.9398306]]
[[1]]
```

Python (Flask):

```
🚵 *app.py - E:\project\Flask\app.py (3.9.6)*
File Edit Format Run Options Window Help
from future import division, print function
import sys
import os
import glob
import numpy as np
from keras.preprocessing import image
from keras.applications.imagenet utils import preprocess input, decode predictio
from keras.models import load model
from keras import backend
from tensorflow.keras import backend
import tensorflow as tf
#global graph
#graph=tf.get default graph()
from skimage.transform import resize
# Flask utils
from flask import Flask, redirect, url_for, request, render_template
from werkzeug.utils import secure filename
from gevent.pywsgi import WSGIServer
import cv2
from tensorflow.python.keras.models import load_model
app = Flask( name )
# Load your trained model
@app.route('/', methods=['GET'])
def index():
    # Main page
    return render template('base.html')
#from skimage.transform import resize
def detect(img path, model):
        #print(type(frame))
        #try:
        img = image.load img(img path, target size = (150, 150))
        img = image.img to array(img)
        img = np.expand dims(img,axis=0)
        if (np.max(img)>1):
            img = img/255.0
        prediction = model.predict(img)
        print (prediction)
                                                                            Ln: 21 Col: 0
```

```
済 *app.py - E:\project\Flask\app.py (3.9.6)*
File Edit Format Run Options Window
                                  <u>H</u>elp
def detect(img path, model):
        #print(type(frame))
        #try:
        img = image.load img(img path, target size = (150, 150))
        img = image.img to array(img)
        img = np.expand dims(img,axis=0)
        if (np.max(img)>1):
            img = img/255.0
        prediction = model.predict(img)
        print (prediction)
        prediction = model.predict classes(img)
        print (prediction)
        return prediction
@app.route('/predict', methods=['GET', 'POST'])
def upload():
    if request.method == 'POST':
        # Get the file from post request
        f = request.files['image']
        # Save the file to ./uploads
        basepath = os.path.dirname( file )
        file path = os.path.join(
            basepath, 'uploads', secure filename(f.filename))
        f.save(file path)
        model = load model("mymodel.h5")
        print('Model loaded. Check http://127.0.0.1:5000/')
        model.compile(loss='categorical crossentropy', optimizer='adam', metrics
        #img = image.load img(file path, target size=(150, 150))
        preds = detect(file path, model)
        #print(float(preds))
        if preds == 0:
            text = "You are perfectly fine"
        else:
            text = "You are infected! Please Consult Doctor"
        text = text
            # ImageNet Decode
        return text
if name == ' main ':
    app.run(debug=True, threaded = False)
                                                                            Ln: 21 Col: 0
```

12.2. UI Output Screenshot

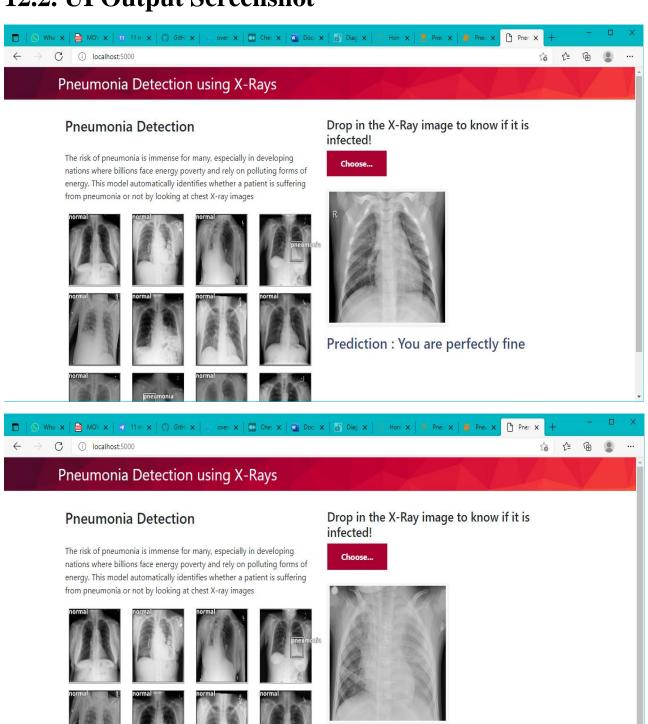

Prediction: You are infected! Please

Consult Doctor